


QTA

Quick Check

- 1 Turn the Tracer (blood bag with Tracer) in your hand until the Tracer light starts flashing.



- 2 Green light - OK.
Red Light - not OK,
blood invalid.

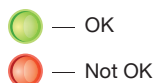

Blood Status and Expiration date

- 1 Start QTA Access Point, click on the "Status" tab.



- 2 Activate QTA Tracer by turning it over in your hand until the light starts flashing.



- 3 Scan the bar code on the Tracer - connection and transfer progress is displayed on screen.



- 4 Current status and expiration date is displayed.



QTA Tracer with Bluetooth 4.0-technology



QTA Tracer System is user-friendly. Here we guide you through activation and handling of your QTA Tracers for best result.

Getting started

- 1 Start QTA Access Point software.



- 2 Scan the bar codes on the blood bag.



- 3 Activate the QTA Tracer by turning it over in your hand until the light starts flashing.



- 4 Scan the bar code on the QTA Tracer. Connection and completion is displayed on screen.



- 5 Attach the QTA Tracer on the blood bag.




Release blood bags from quarantine

- 1 Start QTA Access Point, click on the "Status" tab.



- 2 Activate the QTA Tracer Tracer by turning it over in your hand until the light starts flashing.

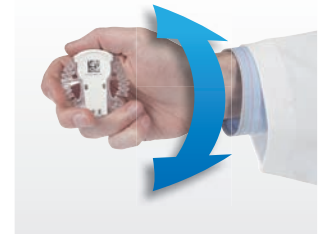

- 3 Scan the bar code on the QTA Tracer. Connection is displayed on screen.



- 4 Click the "release" button or scan bar code on screen to release blood.



- 5 Green light on the QTA Tracer and message on screen confirms the release.

